
Dale® Tracheostomy Tube Holders

APPLICATION INSTRUCTIONS

APPLICATION INSTRUCTIONS

1. Thread long, narrow fastener tabs through flanges on trach plate.
2. Bring tabs back over flange and adhere to the neckband.
3. Secure fastener tab of short section **(A)** to the long neckband section **(B)**. Adjust for proper length or compression.
4. Trim excess material **(B only)** for proper length.

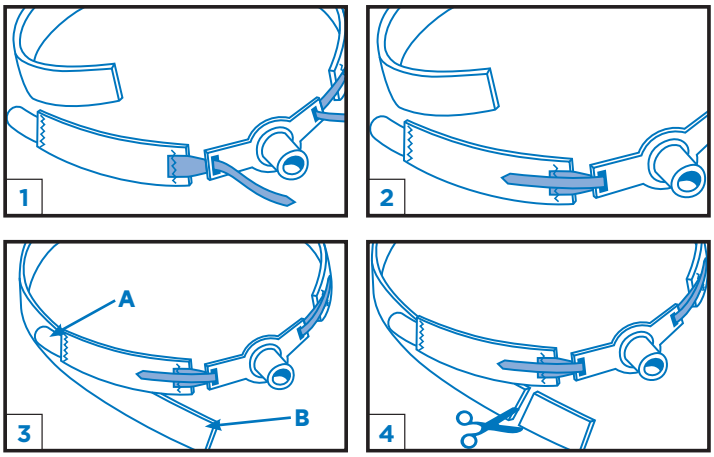

CAUTION

Use in supervised care settings only.

Monitor to prevent unintentional disengagement. Do not trim fastener hook tabs.

RECOMMENDATIONS

Change the holder daily or more frequently as required. Change after patient bathing.

ADDITIONAL BENEFITS

- Secure fastener tabs limit movement of the tracheostomy tube reducing accidental decannulation.
- Exclusive moisture wicking neckband lining reduces the risk of skin breakdown by minimizing the amount of moisture absorbed.
- Adjustable two piece neckband.
- PediPrints® series offers kids fun, comforting prints.†
- Fastener tabs quickly and easily attach to any size tracheostomy tube.
- Stretch material‡ allows for cough reflex, accommodates edema, and ensures a snug fit without compromising security.
- Not made with natural rubber latex.

† On the 241 and 242 tube holder
‡ On the 240 Blue® tube holder



Product #	Description	Quantity
240	**240 Blue®** 2.54cm(1˝) wide, fits up to 49.5cm(19.5˝) neck	10/Box
241	**241 PediPrints®** 2.54cm(1˝) wide, fits up to 45.5cm(18˝) neck	10/Box
242	**242 PediPrints®** 1.91cm(0.75") wide, fits up to 23cm(9˝) neck	10/Box
244	**244 Blue®** 3.81cm(1.5") wide, fits up to 71cm(28˝) neck	10/Box

HCPCS A7526

SCAN FOR MORE INFO

Dale® Medical Products
EMPLOYEE OWNED 100% AND OPERATED

800-343-3980 • www.dalemed.com

Dale, 240 Blue, 244 Blue, & PediPrint are registered trademarks of Dale Medical Products, Inc. in the USA, EU, and UK.
©2024 Dale Medical Products, Inc. All rights reserved.