Dale® BreezeLock®
Endotracheal Tube Holder
APPLICATION INSTRUCTIONS

APPLICATION INSTRUCTIONS

1. Clean tube. Let dry.
2. Apply Adhesive Base.
   - Apply barrier wipe to clean skin & allow to dry.
   - Remove liner and place Base at desired location **(A1-A3)**. Trim if needed.
3. Press neckband onto hook portion of the Adhesive Base, making sure the TRACK is centered above the patients mouth. **(B)**
4. Open the fastener tabs located at both ends of the neckband and wrap the band around the head just below the ears. Secure the fastener tabs to the outside of the neckband.
5. Lift OVERLAY. Position the CHANNEL over the ET Tube. Press OVERLAY across TRACK and Neckband to secure into position. Ensure good contact is made with the neckband, the blue tube track is completely covered, and the plastic channel is vertically centered.
6. Steps to tape the tube to the CHANNEL:
   - with the tube resting freely below the CHANNEL, wrap 8" long piece of strong adhesive tape 1x around the tube. **(C)**
   - Press TUBING CHANNEL down onto the tube.
   - Firmly wrap tape in a forward spiral motion around tube and CHANNEL at least 1" away from patient **(D)**. Be sure to wrap the tube at least 3 times.
7. Ensure neckband security by readjusting the fastener tabs as necessary.

- Change tape on CHANNEL as necessary.
- Change device after 72 hours or more often if needed.
- To prevent skin ulcerations or skin tears, assess patient skin and medical conditions prior to applying. Monitor daily or more frequently.

TUBE REPOSITIONING


While holding the taped CHANNEL with one hand, lift OVERLAY off TRACK. Slide CHANNEL to move tube to desired location; be sure to keep blue track vertically centered. Press OVERLAY back down onto TRACK and neckband to secure the tube's position in place. Ensure good contact is made with the neckband, the blue tube track is completely covered, and the plastic channel is vertically centered.


REMOVAL INSTRUCTIONS

Device should be removed with care to avoid skin trauma. Re-apply skin prep if another device is to be applied.

RECOMMENDATIONS

- **Dale recommends repositioning the endotracheal tube often** to help prevent injury to the lips and underlying tissues due to unrelieved pressure.
- **Dale strongly recommends supporting the ventilator tubing** to reduce pressure on the patient's face.

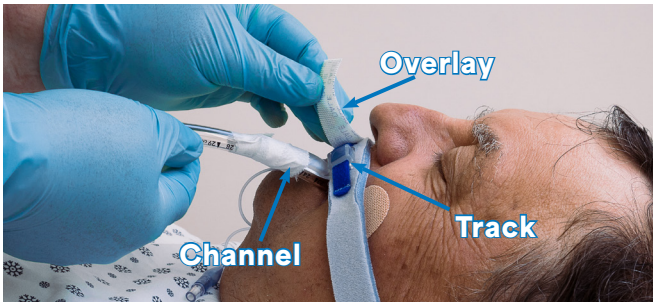

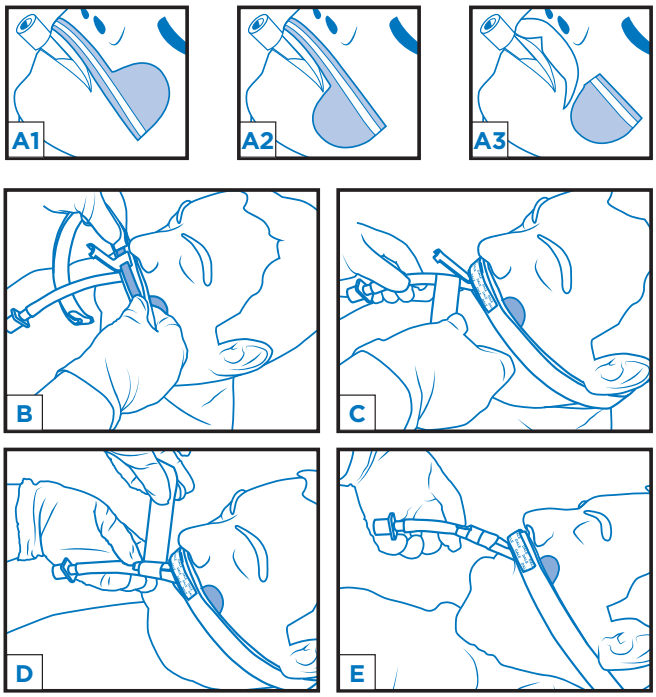

Product #	Description	Quantity
280	**Dale BreezeLock with Tube Track** Fits endotracheal tubes 7.0-10.0mm	10/Box

Visit dalemed.com for an in-service application video

SCAN FOR MORE INFO

Dale Medical Products

EMPLOYEE OWNED • 100% • AND OPERATED

800-343-3980 • www.dalemed.com

Dale & BreezeLock are registered trademarks of Dale Medical Products, Inc. in the USA, EU, and UK.
©2025 Dale Medical Products, Inc. All rights reserved.