REF 280

Dale® BreezeLock® Endotracheal Tube Holder



For Comfortable & Secure Stabilization

Dale® Medical Products

Dale Medical Products, Inc.
40 Kenwood Circle • Suite 7
Franklin, MA 02038-3298, USA
T: +1-508-695-9316
www.dalemed.com

EC REP **MDSS GmbH**
Schiffgraben 41
D-30175 Hannover, Germany

CH REP **MDSS CH GmbH**
Laurenzenvorstadt 61
5000 Aarau, Switzerland

(Croatian, Danish, Dutch, English, Finnish, French, German, Greek, Italian, Japanese, Korean, Norwegian, Polish, Portuguese, Spanish, Swedish, Thai, Turkish)

Dale® Medical Products

BreezeLock® Endotracheal Tube Holder

fits endotracheal tube sizes 7.0-10.0mm
Contents: 1 of each - neckband; 273 Adhesive Base

Application Instructions

1. Clean tube. Let dry.
2. Apply Adhesive Base. **(A)** Follow directions on 273 package.
3. Press neckband onto hook portion of the Adhesive Base, making sure the TRACK is centered above the patients mouth. **(B)**
4. Open the fastener tabs located at both ends of the neckband and wrap the band around the head just below the ears. Secure the fastener tabs to the outside of the neckband.
5. Lift OVERLAY. Position the CHANNEL over the ET Tube. Press OVERLAY across TRACK and Neckband to secure into position. Ensure good contact is made with the neckband, the blue tube track is completely covered, and the plastic channel is vertically centered.
6. Steps to tape the tube to the CHANNEL:
   - With the tube resting freely below the CHANNEL, wrap 8" long piece of strong adhesive tape 1x around the tube. **(C)**
   - Press TUBING CHANNEL down onto the tube.
   - Firmly wrap tape in a forward spiral motion around tube and CHANNEL at least 1" away from patient **(D)**. Be sure to wrap the tube at least 3 times.
7. Ensure neckband security by readjusting the fastener tabs as necessary.

- Change tape on CHANNEL as necessary.
- Change device after 72 hours or more often if needed.
- To prevent skin ulcerations or skin tears, assess patient skin and medical conditions prior to applying. Monitor daily or more frequently.

Tube Repositioning

While holding the taped CHANNEL with one hand, lift OVERLAY off TRACK. Slide CHANNEL to move tube to desired location. Press OVERLAY back down onto TRACK and neckband to secure the tube's position in place. Ensure good contact is made with the neckband, the blue tube track is completely covered, and the plastic channel is vertically centered.

Removal

Device should be removed with care to avoid skin trauma. Re-apply skin prep if another device is to be applied.

Dale **recommends repositioning the endotracheal tube often** to help prevent injury to the lips and underlying tissues due to unrelieved pressure.

Dale **strongly recommends supporting the ventilator tubing** to reduce pressure on the patients face.

For single patient use. Disposable.

For an application video please visit www.dalemed.com


Not made with natural rubber latex.


Patent www.dalemed.com/patents


Notice to the user and/or patient; any serious incident that has occurred in relation to this device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.


Dale & BreezeLock are registered trademarks of Dale Medical Products, Inc. in the USA, EU, and UK.


PHT NON-DEHP DEHP

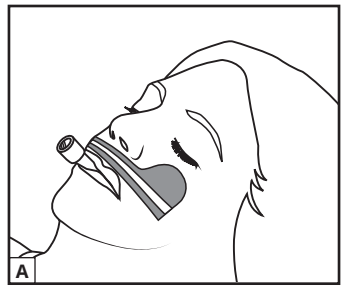

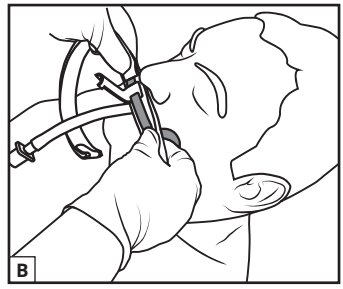

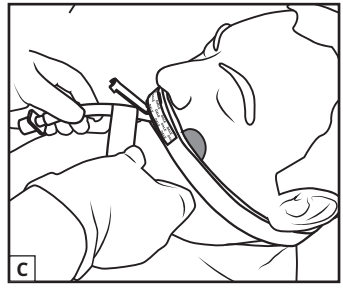

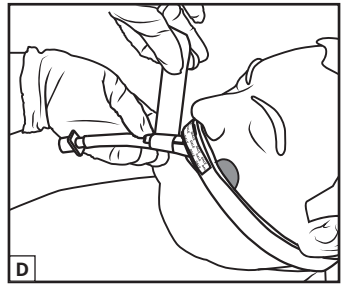

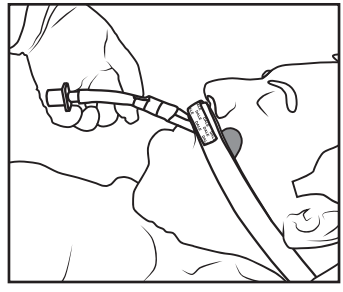

©Copyright 2025

942255 REV C

©Copyright 2025