Dale Hold-n-Place®
Leg Band/ Waist Band
APPLICATION INSTRUCTIONS

316 | LEG BAND APPLICATION INSTRUCTIONS


1. Position leg band under and high around patient's thigh with Dale product label pointing to inside of leg. Stretch leg band equally from both sides into place and fasten hook tab. Proper compression has been achieved when two fingers fit snugly under the band.**(1-2B)**
2. Place Foley catheter over white stitching of green locking device **(3)**. Leave an ample loop in catheter above leg band to avoid bladder traction; or stretch catheter for traction when required.
3. With catheter in desired position, insert the narrow green tab over catheter and through the square opening so the green tab overlaps.**(4)**
4. Lock catheter in position by pulling green tabs in opposite directions. **(5A)** Secure tabs in place. To reposition catheter: raise either green tab, adjust catheter, and refasten tab.**(5B)**
5. Another method of securement is to locate the Y port of the proximal end of the catheter into the locking device.**(6)**


330 | WAIST BAND APPLICATION INSTRUCTIONS


1. To secure the catheter in place, position the waist band around the abdomen, above the wound site. No shaving necessary. The green catheter locking tab should be centered over the drain site. After the band has been positioned, press locking tab onto the waist band.
2. Place catheter over white stitching of green tab.
3. Repeat steps 3,4, & 5 above.


ADDITIONAL BENEFITS


- No shaving required.
- To help prevent meatal irritation, leg band can be rotated in either direction. Alternating legs periodically may also add to greater patient comfort.
- 330 can also be used to secure PEG tubes.

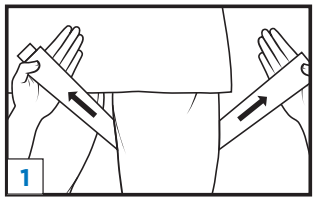


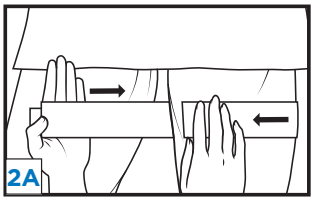


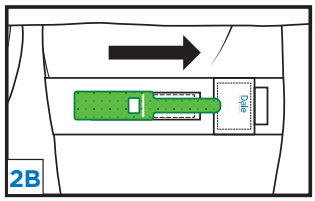


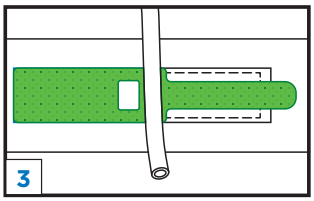

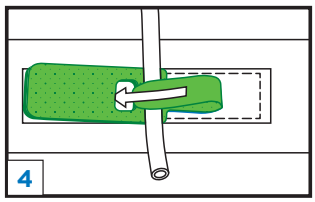


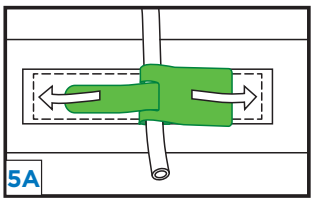

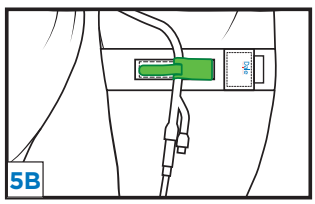

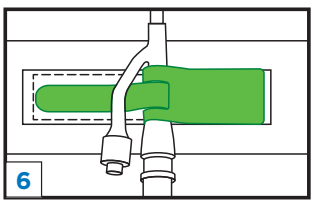



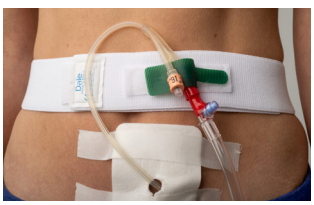

Product #	Description	Quantity
316	Leg band, Fits up to 30"	10/box
330	Leg band/Waist band, Fits up to 56"	10/box

SCAN FOR MORE INFO

Dale Medical Products

EMPLOYEE OWNED 100% AND OPERATED

800-343-3980 • www.dalemed.com

Dale & Hold-n-Place are registered trademarks of Dale Medical Products, Inc.
©2024 Dale Medical Products, Inc. All rights reserved.