
Dale®
Medical Products


Hold-n-Place®
Foley Catheter Tube Holder


Additional Indications - Supra-pubic Foley Catheters or PEG tubes


Waist Band Application


1. Position the waistband around the abdomen, just above the wound site—no shaving required. Ensure the green catheter locking tab is centered over the stoma site. Once the band is in place, firmly press the locking tab onto the waistband to secure it.
2. Place the Y-port of the catheter over the white stitching of the green tab.
3. With catheter in desired position, insert the narrow green tab over catheter and through the square opening so the green tab overlaps.
4. Lock catheter in position by pulling green tabs in opposite directions. Secure tabs in place. To reposition catheter: raise either green tab, adjust catheter, and refasten tab.

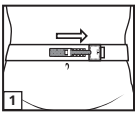

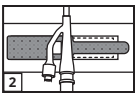

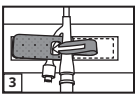

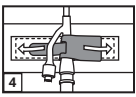

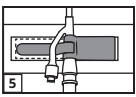

For single patient use. Disposable

For an application video please visit www.dalemed.com

Dale medical products are available from your preferred hospital/medical supply dealer.

Not made with natural rubber latex.

Notice to the user and/or patient; any serious incident that has occurred in relation to this device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

Dale & Hold-n-Place are registered trademarks of Dale Medical Products, Inc. in the USA, EU, and UK

©Copyright 2025