# Dale® ACE Connector® with ENFit® Connection

**Infuse nutrients, medications, fluids, irrigate, check gastric residuals, control flow and suction**

## APPLICATION INSTRUCTIONS

| Product Number |
| --- |
| 485 |


**60mL & most other ENFit sizes**

### ON/OFF HANDLE

controls fluid flow.

"ON" for feeding or drainage.

"OFF" to irrigate tube, check gastric residuals, administer medication, or ambulate patient.

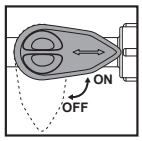


### AUTO-SEAL SYRINGE MED PORT

accepts most 60mL ENFit type tip syringes for irrigation, medication delivery and checking gastric residuals.**

**Align, Push & Twist ENFit 60mL syringe or Administration Set ENFit female connection into ENFit Auto-Seal Syringe Med Port with male ENFit Connection.**

**DO NOT OVERTIGHTEN**

**Medication must be crushed properly to prevent potential clogging.

### PROXIMAL END

**(Feeding Bag End)**

furthest from patient.




### ENFIT MALE CONNECTOR CONNECTS TO:

- enteral pump and gravity delivery sets with ENFit female connection.
- bolus feeding syringes with ENFit tip.
- Suction tube with special ENFit adapter plug.

### DISTAL END

**(Patient End)**

closest to patient's stomach.

- Connects to all enteral tubes with ENFit connectors.

### ENFIT FEMALE CONNECTOR CONNECTS TO:

- Small bore feeding tubes with male ENFit connectors.
- PEG, PEJ, and most G-J tubes with male ENFit connectors.
- Balloon replacement tubes with male ENFit connectors.
- Salem Sump™ tubes equipped with an ENFit connector.

### THE ACE ADVANTAGE

- Supports a closed enteral system
- Protects clinician from potentially infectious splash-back
- Protects patient from outside contamination
- No soiled bed linens or loss of formula due to open adapters
- Provides auto-seal syringe med port seal for accepting medication, irrigation, or gastric residuals

MR CONDITIONAL

Dale Medical Products

EMPLOYEE OWNED • 100% • AND OPERATED

800-343-3980 • www.dalemed.com

SCAN FOR MORE INFO

Note: the distal end is always on the side closest to the pt. This is true from a clinician's perspective. It is exactly opposite from a patient's point of view.

Dale & ACE Connector are registered trademarks of Dale Medical Products, Inc. in the USA, EU, and UK.
ENFit is a registered trademark of GEDSA used with their permission.
Salem Sump is a trademark of Cardinal Health, Inc.
©2024 Dale Medical Products, Inc. All rights reserved.