
Dale® Nasal Dressing Holder
APPLICATION INSTRUCTIONS


APPLICATION INSTRUCTIONS


1. Center plastic section of the holder over the nasal dressing and under the nose. Secure by placing ear loops around patient's ears.
2. To adjust size (*ear to nose*)
   a. Lift end of strap and reposition ear loop assembly to the desired length. Repeat for other side.
   b. Fold end of the strap over the ear loop assembly. Excess strap material can be removed with scissors.

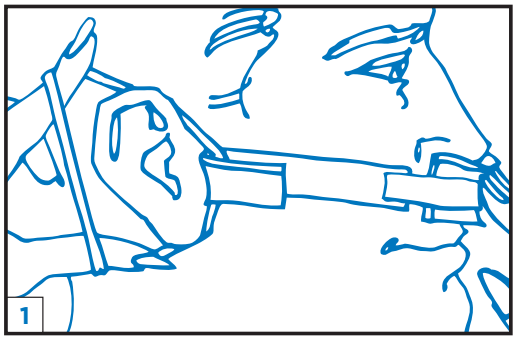

TO CHANGE DRESSING

To change nasal dressing, remove one ear loop, change dressing, and replace ear loop around patient's ear.

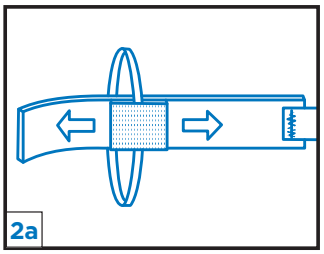

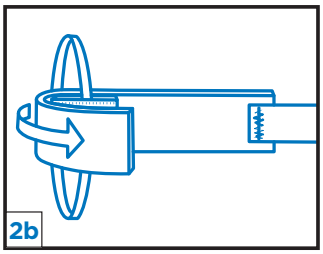

PRODUCT FEATURES

- Eliminates irritation and excoriation of facial skin
- Simplifies dressing changes
- Saves nursing time
- Improves patient comfort
- One size fits most patients
- Not made with natural rubber latex



Product #	Description	Quantity
600	One size fits most	10/Box

SCAN FOR MORE INFO

Dale Medical Products

800-343-3980 • www.dalemed.com

Dale is a registered trademark of Dale Medical Products, Inc. in the USA, EU, and UK.
©2024 Dale Medical Products, Inc. All rights reserved.