

Dale® Transducer Holder
A Level and Secure Platform
APPLICATION INSTRUCTIONS


APPLICATION INSTRUCTIONS

1. Position the holder along the phlebostatic axis (*mid-chest*) and wrap around the patient's arm. Press hook & loop closure to secure.
2. Open the necessary hook & loop tabs, thread through the slot(*s*) of the transducer(*s*) and press to secure.



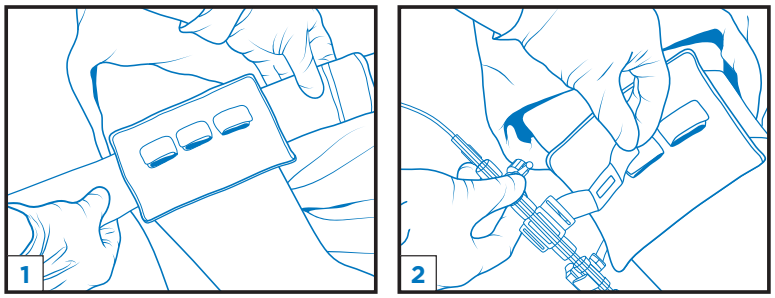

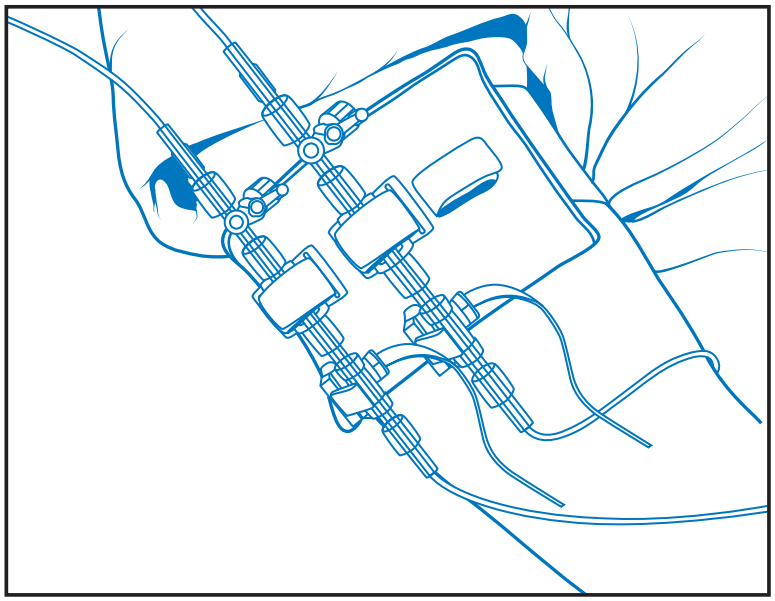

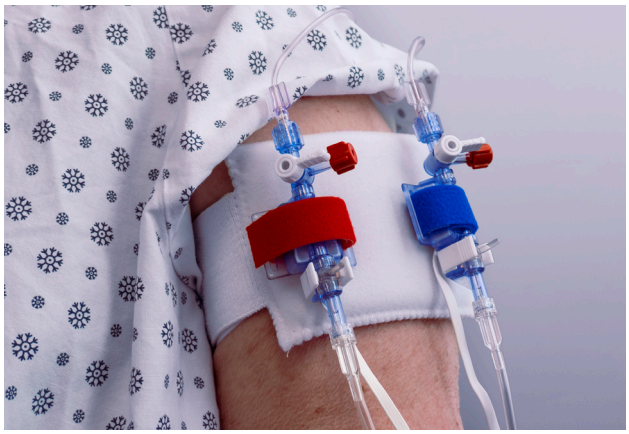

Product #	Description	Quantity
620	**Fits up to three standard transducers** One size fits most	10/Box
622	**Fits up to two large transducers** One size fits most	10/Box

PRODUCT FEATURES

- Soft platform secures and organizes up to three transducers.
- Allows for a complete securement system for radial artery blood pressure monitoring when used in conjunction with the Dale Bendable ArmBoard.
- Quick and easy to apply and remove. Can be used over BP Cuff or sleeve.
- Saves nursing time by eliminating the need to tape individual transducers to the patient's chest, arm or mounting on bedside pole.
- One size fits most.
- Not made with natural rubber latex.

SCAN FOR MORE INFO

Dale® Medical Products

EMPLOYEE OWNED • 100% • AND OPERATED

800-343-3980 • www.dalemed.com

Dale is a registered trademark of Dale Medical Products, Inc. in the USA, EU, and UK.
©2024 Dale Medical Products, Inc. All rights reserved.