
Dale® IV-ARMOR®
APPLICATION INSTRUCTIONS

APPLICATION INSTRUCTIONS


1. Secure and dress IV line per hospital protocol.
2. Apply skin prep evenly to areas adhesive will adhere to, and let dry for approximately 30 seconds.
3. Center the pad over the IV & dressing. Determine if length of legs need to be cut to prevent overlapping. If necessary, cut legs before removing liner.

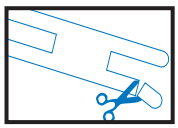



*Do not cut more than 2" from one side. If you require cutting more than 2", you must trim from other set of product legs.

4. Pull back liner and press product to adhere to patient in one motion, one side at a time. Do not stretch to apply.




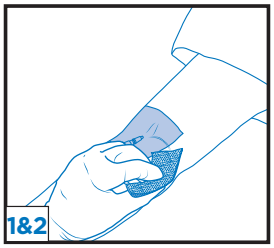

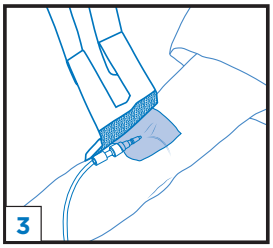


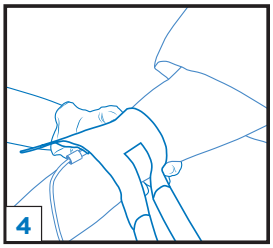


ANTECUBITAL APPLICATION

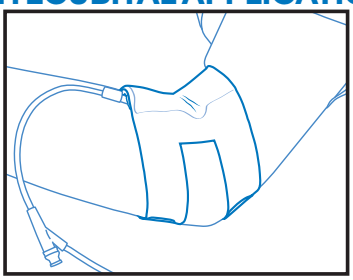

HAND APPLICATION

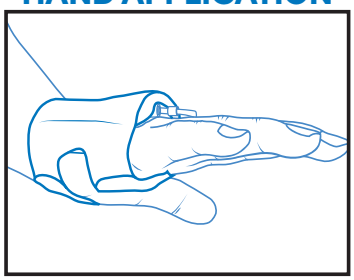

FOREARM APPLICATION

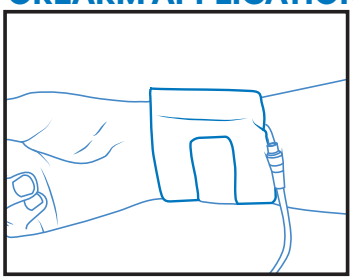

Keep device clean and dry.
Change after 3 days or sooner if necessary.

PRODUCT FEATURES

- Helps maintain a continuous IV flow while protecting and preserving the IV site
- Reduces the number of IV restarts, which decreases risk of infection
- Saves nursing time and reduces alarm fatigue
- Allows for patient movement while reducing downstream occlusion alarms caused by movement
- Improves patient comfort without the need for a stiff immobilizing device
- Reduces risk of dislodgement from picking and pulling of the IV by confused patients

REMOVAL INSTRUCTIONS

Device should be removed with care to avoid skin trauma. Re-apply skin prep if another device is to be applied.

Product #	Description	Quantity
675	One size fits most	10/Box

HCPCS Codes A4221

SCAN FOR MORE INFO

Dale Medical Products
EMPLOYEE OWNED • 100% • AND OPERATED

800-343-3980 • www.dalemed.com

Dale & IV-ARMOR are registered trademarks of Dale Medical Products, Inc. in the USA, EU, and UK.
©2024 Dale Medical Products, Inc. All rights reserved.