
Dale® Post-Surgical Bra
APPLICATION INSTRUCTIONS

APPLICATION INSTRUCTIONS

1. Detach straps and place the bra under the patient.
2. Grasp both sides of the bra and press together. There should be some tension.
3. Attach straps to front of bra at the top of the band.
4. Bra can easily be adjusted by peeling fasteners apart at front and straps.

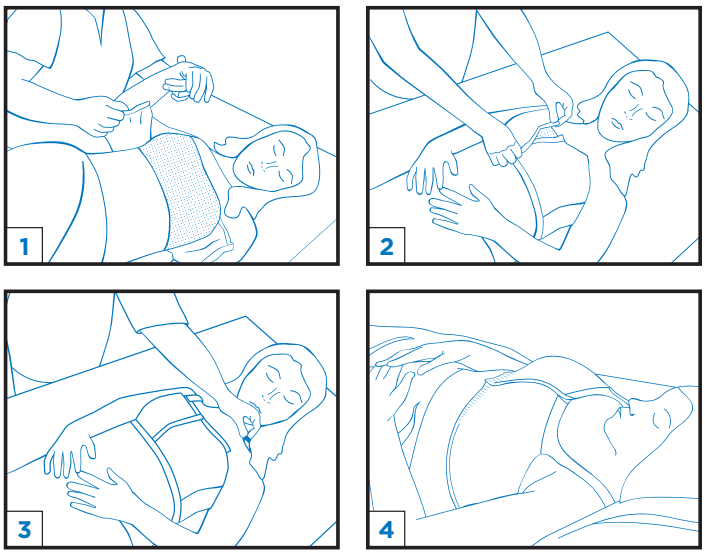



ALTERNATIVE APPLICATION

The Dale Post Surgical Bra may also be applied without detaching the straps by gently guiding the patient's arms through the straps.

PRODUCT FEATURES

- Secure, long-term comfort
- Stretch knit fabric lets treatment site breathe
- Material expands and contracts to accommodate edema
- Seamless design helps prevent trauma
- Hook and loop closure permits easy application and inspection
- Detachable straps allow for multiple application options
- Provides maximum comfort throughout recovery
- Not made with natural rubber latex.

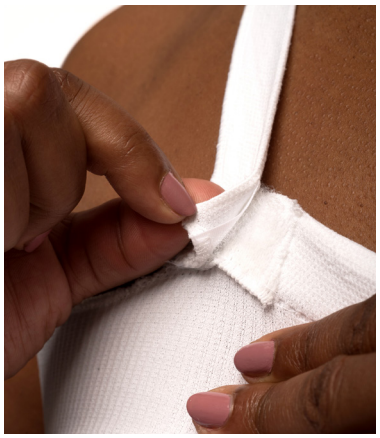

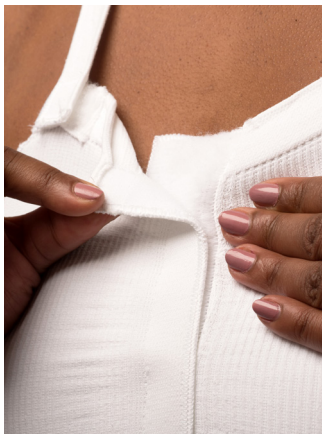



Product #	Description	Quantity
701	Small, fits A-C 81cm-86cm (32"-34")	1/Box
702	Medium, fits B-D 86cm-92cm (34"-36")	1/Box
703	Large, fits B-D 92cm-97cm (36"-38")	1/Box
704	X-Large, fits B-D 97cm-112cm (38"-44")	1/Box
705	XX-Large, fits C-E 117cm-137cm (46"-54")	1/Box
706	XXX-Large, fits E-I 137cm-152cm (54"-60")	1/Box

HCPCS Codes A4467

SCAN FOR MORE INFO

Dale Medical Products

EMPLOYEE OWNED • 100% • AND OPERATED

800-343-3980 • www.dalemed.com

Dale is a registered trademark of Dale Medical Products, Inc. in the USA, EU, and UK.
©2026 Dale Medical Products, Inc. All rights reserved.